

Reflectance Spectra of Aluminum Mirror

Sample


Sample : Aluminum Flat Mirror TFAN-30S05-1
(Sigma Koki, Tokyo, Japan)

Instrument condition

INSTRUMENT	: U-4100		
UV/VIS		NIR	
SCAN SPEED	: 300 nm/min	SCAN SPEED	: 750 nm/min
BANDPASS	: 8 nm	BANDPASS	: AUTO CHANGE
		PbS Gain	: 2

Accessory

45° specular reflectance accessory
(P/N : 134-0106)
Polarizer holder (P/N : 132-0325)
Glan-Taylor polarizing prism
MGTYB20(Karl Lambrecht, Chicago,
USA)

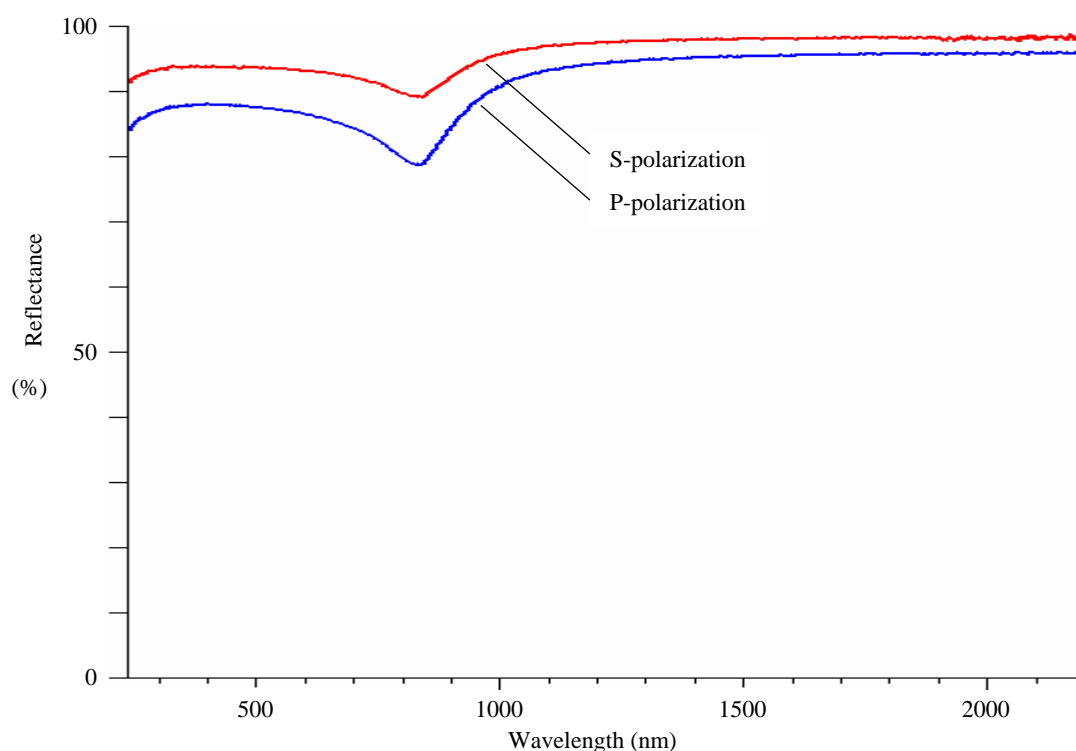

KEY WORDS

Aluminum Mirror, Reflectance, 45° Specular Reflectance, Reflectance Spectrum, UV, U-4100

Spectrophotometer UV

Sheet No. UV070006-01